

# **Spectroradiometer CS-2000 / CS-2000A**

## **Communication Specifications**



**KONICA MINOLTA**

Konica Minolta, Inc.  
Optics Company, Sensing Business Unit  
3-91, Daisennishimachi, Sakai-ku, Sakai-shi, Osaka, JAPAN  
<http://www.konicaminolta.com/selector/instruments.html>

© 2009-2014 KONICA MINOLTA, INC.

**Notes regarding these specifications:**

- The commands and other information contained within this document cannot be used in commercial software without the express permission of Konica Minolta, Inc.
- An understanding of the fundamentals of PC communication is assumed in the preparation of this document. This document was prepared to explain communication procedures specific to the stated product.
- Reprinting of this document or the information contained herein in its entirety or in part without the express permission of Konica Minolta, Inc. is strictly prohibited.
- The information herein is subject to change without notice.
- Although every effort has been made to ensure the accuracy of the information contained herein, if you discover an error, missing information, etc. please notify your Konica Minolta representative.
- Konica Minolta is not responsible in any way for any losses, damage, etc. resulting from the use of the information contained herein.
- Company names and product names contained herein are trademarks or registered trademarks of their respective companies.

Revision (Date)	Modifications
Rev. 1.00 (March 13, 2009)	First English version
Rev. 1.01 (September 12, 2013)	Change of company name
Rev. 1.02 (December 8, 2014)	Minor text edits

## Table of Contents

<b>1</b>	<b>Communication Overview .....</b>	<b>5</b>
1.1	Communication method.....	5
1.1.1	USB.....	5
1.2	Instrument control.....	5
1.3	CS-2000/CS-2000A settings (Remote Mode).....	5
1.3.1	Operation of instrument keys in Remote Mode.....	6
1.3.2	Enabling instrument measurement button.....	6
1.4	Delimiter .....	6
1.5	Data string formats .....	6
1.6	Timeout setting.....	7
1.7	Storage of settings.....	7
1.8	Hexadecimal format used in commands.....	7
<b>2</b>	<b>Commands.....</b>	<b>8</b>
2.1	Command list.....	8
<b>3</b>	<b>Command details .....</b>	<b>10</b>
	RMTS (Remote Mode Select).....	10
	IDDR (Identification Data Read) .....	11
	SCMR (Sync Mode Read) .....	12
	SCMS (Sync Mode Set).....	13
	SPMR (Speed Mode Read) .....	14
	SPMS (Speed Mode Set).....	16
	STSR (Aperture Stop Status Read) .....	18
	UCCS (User Calibration Channel Select) .....	19
	UCCR (User Calibration Channel Read) .....	20
	UCPS (User Calibration Parameter Set).....	21
	UCPR (User Calibration Parameter Read) .....	23
	UCCD (User Calibration Channel Delete).....	24
	LNSS (Close-up Lens Status Select).....	25
	LNSR (Close-up Lens Status Read) .....	26
	ALFS (Attachment Lens Compensation Factor Set).....	27
	ALFR (Attachment Lens Compensation Factor Read) .....	28
	NDFS (External ND Filter Select) .....	29
	NDFR (External ND Filter Read).....	30
	NFCS (ND Filter Compensation Factor Set).....	31
	NFCR (ND Filter Compensation Factor Read) .....	32
	MEAS (Measure) .....	33
	MEDR (Measurement Data Read).....	35
	MSWE (Measuring Switch Enable).....	38
	STDS (Store Data Set) .....	39
	STDR (Stored Measurement Data Read) .....	40
	STDD (Store Data Delete) .....	42
	STAD (Store All Data Delete) .....	43
	TGSL (Target Number Select).....	44
	TGSR (Target Number Selection Read) .....	45
	TGDS (Target Data Set) .....	46
	TGDR (Target Data Read).....	48
	TGDD (Target Data Delete).....	51
	TGAD (Target All Data Delete) .....	52

BALS (Backlight Control Set).....	53
BALR (Backlight Control Read) .....	54
CSMS (Color Space Mode Set) .....	55
CSMR (Color Space Mode Read).....	56
DIMS (Display Mode Set) .....	57
DIMR (Display Mode Read) .....	58
OBSS (Observer Set) .....	59
OBSR (Observer Read) .....	60
DTCR (Date/Time of Calibration Read) .....	61
<b>4 Error-check codes.....</b>	<b>62</b>
<b>5 Numerical output formats.....</b>	<b>63</b>
5.1 Spectral irradiance values .....	63
5.1.1 Alphanumeric data .....	63
5.1.2 Hexadecimal data.....	63
5.2 Colorimetric values .....	64
5.2.1 Alphanumeric data .....	65
5.2.2 Hexadecimal data.....	67
<b>6 Numerical input formats .....</b>	<b>68</b>
6.1 Spectral irradiance values .....	68
6.1.1 Alphanumeric data .....	68
6.1.2 Hexadecimal data.....	68
6.2 Colorimetric values .....	69
6.2.1 Alphanumeric data .....	70
6.2.2 Hexadecimal data.....	70
<b>7 Characters for data communication.....</b>	<b>71</b>
<b>8 Measurement flow .....</b>	<b>72</b>
8.1 Communication flow for measurement controlled from PC.....	72
8.2 Communication flow for measurement using instrument measuring key..	73
<b>9 Modification of programs written for CS-2000 units with previous firmware versions.....</b>	<b>74</b>
9.1 Function changes in updated firmware .....	74
9.2 Command changes to correspond to updated firmware .....	74
9.2.1 Command changes to correspond to change of output for calculation error number .....	74
9.2.2 Command changes to correspond to calculation error number changes.....	75
9.2.3 Command changes to correspond to Identification Data changes.....	75

# 1 Communication Overview

This communication protocol is for both Konica Minolta Spectroradiometer CS-2000 and CS-2000A.

## 1.1 Communication method

Communication with the CS-2000/CS-2000A is performed via USB.

- For information on connecting the CS-2000/CS-2000A to a computer, please refer to the “Communications” section of the Spectroradiometer CS-2000/CS-2000A Instruction Manual.

### 1.1.1 USB

Communication with a PC can be performed with the CS-2000/CS-2000A configured as using a COM port. The settings are shown below. Note that communication can be performed regardless of the baud rate setting.

Specification		USB 1.1 Full Speed
File used		Kmse2000.inf (Installed into CS-S10w program folder when CS-S10w software included with the instrument as a standard accessory is installed.)
Communication settings	Baud rate	9,600
	Data length	8 bits
	Parity	None
	Stop bits	1
	Flow control	None

## 1.2 Instrument control

Control of the CS-2000/CS-2000A is performed by sending the command and any associated command parameters for the desired operation from the PC to the CS-2000/CS-2000A.

After the requested operation has been performed, the CS-2000/CS-2000A sends back the operation results as an error-check code and any associated response parameters.

- Command details (command parameters, error-check codes, and response parameters) are explained in the Command Details section.

## 1.3 CS-2000/CS-2000A settings (Remote Mode)

When controlling the CS-2000/CS-2000A from a PC, it is necessary to set the instrument to Remote Mode.

Procedure for setting Remote Mode

- 1 Connect PC and instrument with cable.

- 2 Send Remote Mode setting command "RMTS,1" from PC to instrument.
- 3 Check that error-check code "OK00" is returned from the instrument.

### 1.3.1 Operation of instrument keys in Remote Mode

When the instrument is set to Remote Mode, communication commands can be used to operate the instrument, and instrument keys except for ESC are disabled. In Remote Mode, the instrument ESC key has the following functions:

- When the instrument is not taking a measurement, pressing the instrument ESC key will cancel Remote Mode.
- When the instrument is measuring, pressing the instrument ESC key will interrupt and cancel the measurement.

### 1.3.2 Enabling instrument measurement button

The measurement button of the instrument can be enabled so that it can be used even in Remote Mode using the command "MSWE".

When the measurement button has been enabled in Remote Mode, measurements can be started by either sending the command "MEAS" or by pressing the instrument's measurement button.

When the measurement button has been enabled, the timing for clearing of measurement data stored in the CS-2000/CS-2000A is different. For details, refer to the Explanation section for the command "MEDR".

## 1.4 Delimiter

When sending commands from a PC to the instrument, any of the following delimiters can be used:

[CR] (Carriage return: 0x0D)

[LF] (Line feed: 0x0A)

[CR] (Carriage return: 0x0D) + [LF] (Line feed: 0x0A)

When the instrument returns the error-check code and/or data to the PC in response to a command, the delimiter code used by the instrument will be the same one that was used with the command.

If there is no delimiter code, the command will not be properly recognized. In such case, if nothing is received by the instrument within 60s after receipt of the last character, the received character string will be cleared and the command will be considered invalid.

(Note 1)

For CS-2000 firmware versions of 1.01.0000, if nothing is received by the instrument within 60s after receipt of the last character, the received character string will be cleared and the error-check code "ER00" will be returned.

## 1.5 Data string formats

When sending numerical value, the number of numerical characters in the string must not exceed the number of characters indicated for the command, but can be fewer than the indicated number.

When sending a text string, the number of characters in the string must be the same as the number of characters indicated for the command. Therefore, if the number of characters in the desired string is fewer than the indicated number, spaces should be added to the initial string to create a string of the required length.

## **1.6 Timeout setting**

The timeout time for the PC communication port should be set to at least 10s.

## **1.7 Storage of settings**

The measurement conditions, compensation factors, target values, and instrument display information set using commands are stored in the internal memory of the CS-2000/CS-2000A. As a result, such settings are maintained even after the instrument Remote Mode has been canceled.

## **1.8 Hexadecimal format used in commands**

The parameters of some commands are specified as being in hexadecimal format. The hexadecimal format for use with the CS-2000/CS-2000A is IEEE floating point format (4-byte big-endian hexadecimal string)

## 2 Commands

This communication protocol is for both the CS-2000 and CS-2000A.

However, for the CS-2000, there are some differences in command input and parameter input/output for firmware versions of 1.10.0003 and later compared to versions 1.01.0000 or earlier. For details, refer to the command explanations.

- Parameters shown in parentheses in the “Input/Output Format” should not be input or are not output in some cases, depending on other parameters. Please see the “Command Parameters” and “Response Parameters” sections of each command for details.

### 2.1 Command list

The commands for the CS-2000 and CS-2000A are shown in the table below. Commands for which input/returned parameters for firmware ver. 1.10.0003 are different from those of firmware ver. 1.01.0000 are indicated by “Y” in the Changed column.

Command	Description	Changed	Page
Instrument information/Status			
RMTS	Remote Mode Select		10
IDDR	Identification Data Read	Y	11
Condition settings			
SCMR	Sync Mode Read		12
SCMS	Sync Mode Set		13
SPMR	Speed Mode Read	Y	14
SPMS	Speed Mode Set	Y	16
STSR	Aperture Stop Status Read		18
Calibration settings			
UCCS	User Calibration Channel Select		19
UCCR	User Calibration Channel Read		20
UCPS	User Calibration Parameter Set		21
UCPR	User Calibration Parameter Read		23
UCCD	User Calibration Channel Delete		24
Optional Close-Up Lens settings			
LNSS	Close-up Lens Status Select		25
LNSR	Close-up Lens Status Read		26
ALFS	Attachment Lens Compensation Factor Set		27
ALFR	Attachment Lens Compensation Factor Read		28
Optional ND Filter settings			
NDFS	External ND Filter Select		29
NDFR	External ND Filter Read		30
NFCS	ND Filter Compensation Factor Set		31
NFCR	ND Filter Compensation Factor Read		32
Measurement			
MEAS	Measure		33
MEDR	Measurement Data Read	Y	35
MSWE	Measuring Switch Enable		38
Measurement values			
STDS	Store Data Set		39



STDR	Stored Measurement Data Read	Y	40
STDD	Store Data Delete		42
STAD	Store All Data Delete		43
Target colors			
TGSL	Target Number Select		44
TGSR	Target Number Selection Read		45
TGDS	Target Data Set		46
TGDR	Target Data Read	Y	48
TGDD	Target Data Delete		51
TGAD	Target All Data Delete		52
Settings			
BALS	Backlight Control Set		53
BALR	Backlight Control Read		54
CSMS	Color Space Mode Set		55
CSMR	Color Space Mode Read		56
DIMS	Display Mode Set		57
DIMR	Display Mode Read		58
OBSS	Observer Set		59
OBSR	Observer Read		60

### 3 Command details

RMTS ( <u>R</u> emote <u>M</u> ode <u>S</u> elect)			
Function			
Selects the remote mode setting: On or Off.			
Input/Output Format			
PC		CS-2000/CS-2000A	
"RMTS," + <span style="border: 1px dashed black; padding: 0 5px;">1</span> + <span style="border: 1px dashed black; padding: 0 5px;">Delimiter code</span>		⇒	
		⇐	" <span style="border: 1px dashed black; padding: 0 5px;">Error-check code</span> " + <span style="border: 1px dashed black; padding: 0 5px;">Delimiter code</span>
Command Parameters			
	Meaning	Type	Details/range
<span style="border: 1px dashed black; padding: 0 5px;">1</span>	Remote mode setting	Integer; 1 digit	0: Off (Key mode: commands other than RMTS are not accepted.) 1: On (Communication commands will be accepted.)
Response Parameters			
	Meaning	Type	Details/range
Explanation			
<p>Sets the remote mode setting.</p> <p>When remote mode is off (0), the unit is in key mode and will not accept any commands other than RMTS. (If other commands are sent, the instrument will return "ER00".)</p> <p>When remote mode is on (1), the unit will accept communication commands, and most key operations are disabled, with the following exceptions:</p> <ul style="list-style-type: none"> <li>• When the instrument ESC key is pressed while the instrument is not taking a measurement, remote mode will be canceled.</li> <li>• When the instrument ESC key is pressed while the unit is taking a measurement, the measurement will be canceled.</li> <li>• When the instrument measurement button has been enabled (using the command <u>MSWE</u>), pressing the measurement button starts measurements.</li> </ul>			
Error-check codes			
Code	Meaning		
OK00	Normal completion		
ER00	Invalid command string/number of parameters		
ER17	Invalid parameter received.		

IDDR ( <u>I</u> dentification <u>D</u> ata <u>R</u> ead)			
Function			
Reads the product identification information from the instrument.			
Input/Output Format			
PC		CS-2000/CS-2000A	
"IDDR" + <span style="border: 1px dashed black; padding: 2px;">Delimiter code</span>		⇒	
		⇐	" <span style="border: 1px dashed black; padding: 2px;">Error-check code</span> , <span style="border: 1px dashed black; padding: 2px;">1</span> , <span style="border: 1px dashed black; padding: 2px;">2</span> , <span style="border: 1px dashed black; padding: 2px;">3</span> "
		+	<span style="border: 1px dashed black; padding: 2px;">Delimiter code</span>
Command Parameters			
	Meaning	Type	Details/range
Response Parameters			
	Meaning	Type	Details/range
<span style="border: 1px dashed black; padding: 2px;">1</span>	Product name	String; 9 char.	"CS-2000" or "CS-2000A" May be different for custom units.
<span style="border: 1px dashed black; padding: 2px;">2</span>	Variation code	Integer; 1 digit	1: CS-2000 2: CS-2000A 0 to 9: Used to differentiate models Meaning of numbers other than "1" and "2" is undefined.
<span style="border: 1px dashed black; padding: 2px;">3</span>	Serial number	Integer; 7 digits	
Explanation			
Reads the product identification information (product name, variation, serial number) from the instrument.			
Error-check codes			
Code	Meaning		
OK00	Normal completion		

SCMR (Sync Mode Read)			
Function			
Reads currently set sync (synchronization) mode of the instrument.			
Input/Output Format			
PC		CS-2000/CS-2000A	
"SCMR" + <span style="border: 1px solid black; padding: 2px;">Delimiter code</span>		⇒	
		⇐	" <span style="border: 1px solid black; padding: 2px;">Error-check code</span> , <span style="border: 1px solid black; padding: 2px;">[1]</span> (, <span style="border: 1px solid black; padding: 2px;">[2]</span> )"
		+	<span style="border: 1px solid black; padding: 2px;">Delimiter code</span>
Command Parameters			
	Meaning	Type	Details/range
Response Parameters			
	Meaning	Type	Details/range
<span style="border: 1px solid black; padding: 2px;">[1]</span>	Sync mode	Integer; 1 digit	0: No sync 1: Internal sync 2: External sync
<span style="border: 1px solid black; padding: 2px;">[2]</span>	When Response Parameter [1] = "0" (No sync) or "2" (External sync): Not output		
	When Response Parameter [1] = "1" (Internal sync): Synchronization frequency	Integer; 5 digits	2000 to 20000: Synchronization frequency for internal sync. 100x actual value (Actual range: 20.00 to 200.00Hz) • If number of digits is fewer than 5, "0" will be added before value.
Explanation			
Reads current instrument Sync mode. When Sync mode is set to 1 (Internal sync), a second parameter indicating the synchronization frequency will also be output.			
Error-check codes			
Code	Meaning		
OK00	Normal completion		



SPMR ( <u>Speed</u> <u>Mode</u> <u>Read</u> )			
Function			
Reads currently set speed mode of the instrument.			
Input/Output Format			
For CS-2000 with instrument firmware ver. 1.01.0000 or earlier			
PC		CS-2000	
"SPMR" + <span style="border: 1px dashed black; padding: 2px;">Delimiter code</span>		⇒	
		⇐ " <span style="border: 1px dashed black; padding: 2px;">Error-check code</span> , <span style="border: 1px dashed black; padding: 2px;">[1]</span> ( , <span style="border: 1px dashed black; padding: 2px;">[2]</span> , <span style="border: 1px dashed black; padding: 2px;">[3]</span> ) "	
		+ <span style="border: 1px dashed black; padding: 2px;">Delimiter code</span>	
For CS-2000 with instrument firmware ver. 1.10.0003 or later, or for CS-2000A			
PC		CS-2000/CS-2000A	
"SPMR" + <span style="border: 1px dashed black; padding: 2px;">Delimiter code</span>		⇒	
		⇐ " <span style="border: 1px dashed black; padding: 2px;">Error-check code</span> , <span style="border: 1px dashed black; padding: 2px;">[1]</span> , <span style="border: 1px dashed black; padding: 2px;">[2]</span> ( , <span style="border: 1px dashed black; padding: 2px;">[3]</span> ) "	
		+ <span style="border: 1px dashed black; padding: 2px;">Delimiter code</span>	
Command Parameters			
	Meaning	Type	Details/range
Response Parameters			
For CS-2000 with instrument firmware ver. 1.01.0000 or earlier			
	Meaning	Type	Details/range
<span style="border: 1px dashed black; padding: 2px;">[1]</span>	Speed mode	Integer; 1 digit	0: Normal                      1: Fast 2: Multi Integ                3: Manual
<span style="border: 1px dashed black; padding: 2px;">[2]</span>	<u>When Response Parameter [1] = "0" (Normal) or "1" (Fast):</u> Not output		
	<u>When Response Parameter [1] = "2" (Multi Integ):</u> Integration time	Integer; 2 digits	01 to 16: Integration time in seconds • If number of digits is fewer than 2, "0" will be added before value.
	<u>When Response Parameter [1] = "3" (Manual):</u> Integration time	Integer; 9 digits	000005000 to 120000000: Measurement time in µsec. • If number of digits is fewer than 9, "0" will be added before value.
<span style="border: 1px dashed black; padding: 2px;">[3]</span>	<u>When Response Parameter [1] = "0" (Normal), "1" (Fast), or "2" (Multi Integ):</u> Not output		
	<u>When Response Parameter [1] is set to "3" (Manual):</u> Internal ND filter mode	Integer; 1 digit	0: Off 1: On

For CS-2000 with instrument firmware ver. 1.10.0003 or later, or for CS-2000A			
	Meaning	Type	Details/range
[1]	Speed mode	Integer; 1 digit	0: Normal                      1: Fast 2: Multi Integ Normal      3: Manual 4: Multi Integ Fast
[2]	<u>When Response Parameter [1] is set to "0" (Normal) or "1" (Fast):</u> Internal ND filter mode	Integer; 1 digit	0: Off                              1: On 2: Auto
	<u>When Response Parameter [1] is set to "2"(Multi Integ Normal) or "4"(Multi Integ Fast):</u> Integration time	Integer; 2 digits	01 to 16: Integration time in seconds • If number of digits is fewer than 2, "0" will be added before value.
	<u>When Response Parameter [1] is set to "3" (Manual):</u> Integration time	Integer; 9 digits	000005000 to 120000000: Integration time in $\mu$ sec. • If number of digits is fewer than 9, "0" will be added before value.
[3]	<u>When Response Parameter [1] = "0" (Normal) or "1" (Fast):</u> Not output		
	<u>When Response Parameter [1] is set to "2"(Multi Integ Normal) or "4"(Multi Integ Fast):</u> Internal ND filter mode	Integer; 1 digit	0: Off                              1: On 2: Auto
	<u>When Response Parameter [1] is set to "3" (Manual):</u> Internal ND filter mode	Integer; 1 digit	0: Off                              1: On
<b>Explanation</b>			
Reads Speed mode, Integration time, and Internal ND Filter Mode currently set on instrument. <ul style="list-style-type: none"> <li>For CS-2000 with instrument firmware ver. 1.01.0000 or earlier, Internal ND Filter Mode is read only when Speed Mode is set to Manual (Response Parameter [1] is set to "3"); for other Speed Modes, Internal ND Filter Mode is fixed at Auto.</li> </ul>			
<b>Error-check codes</b>			
Code	Meaning		
OK00	Normal completion		

SPMS ( <u>S</u> peed <u>M</u> ode <u>S</u> et)			
Function			
Sets speed mode of the instrument.			
Input/Output Format			
PC		CS-2000/CS-2000A	
"SPMS, [1] (, [2] [3] )" + [Delimiter]		⇒	
[code]		⇐ "[Error-check code]" + [Delimiter]	
		[code]	
Command Parameters			
For CS-2000 with instrument firmware ver. 1.01.0000 or earlier			
	Meaning	Type	Details/range
[1]	Speed mode	Integer; 1 digit	0: Normal                      1: Fast 2: Multi Integ                3: Manual
[2]	When Command Parameter [1] = "0" (Normal) or "1" (Fast): Do not input.		<ul style="list-style-type: none"><li>If a value is input, "ER00" (Invalid command string or number of parameters) will be returned.</li></ul>
	When Command Parameter [1] is set to "2" (Multi Integ): Integration time	Integer; 2 digits	1 to 16: Integration time in seconds
	When Command Parameter [1] is set to "3" (Manual): Integration time	Integer; 9 digits	5000 to 120000000: Measurement time in μsec.
[3]	When Command Parameter [1] = "0" (Normal), "1" (Fast), or "2" (Multi Integ): Do not input.		<ul style="list-style-type: none"><li>If a value is input, "ER00" (Invalid command string or number of parameters) will be returned.</li></ul>
	When Command Parameter [1] is set to "3" (Manual): Internal ND filter mode	Integer; 1 digit	0: Off 1: On
For CS-2000 with instrument firmware ver. 1.10.0003 or later, or for CS-2000A			
	Meaning	Type	Details/range
[1]	Speed mode	Integer; 1 digit	0: Normal                      1: Fast 2: Multi Integ Normal        3: Manual 4: Multi Integ Fast
[2]	When Command Parameter [1] is set to "0" (Normal) or "1" (Fast): Internal ND filter mode	Integer; 1 digit	0: Off                              1: On 2: Auto <ul style="list-style-type: none"><li>If not set, will automatically be set to "2" (Auto)</li></ul>
	When Command Parameter [1] is set to "2"(Multi Integ Normal) or "4"(Multi Integ Fast): Integration time	Integer; 2 digits	1 to 16: Integration time in seconds
	When Command Parameter [1] is set to "3" (Manual): Integration time	Integer; 9 digits	5000 to 120000000: Measurement time in μsec.



Page 17/75

STSR (Aperture <u>S</u> top <u>S</u> tatus <u>R</u> ead)			
Function			
Reads the status of the instrument's aperture stop (measurement angle).			
Input/Output Format			
PC		CS-2000/CS-2000A	
"STSR" + <span style="border: 1px dashed black; padding: 2px;">Delimiter code</span>		⇒	
		⇐	" <span style="border: 1px dashed black; padding: 2px;">Error-check code</span> , <span style="border: 1px dashed black; padding: 2px;">I</span> " + <span style="border: 1px dashed black; padding: 2px;">Delimiter</span> <span style="border: 1px dashed black; padding: 2px;">code</span>
Command Parameters			
	Meaning	Type	Details/range
Response Parameters			
	Meaning	Type	Details/range
<span style="border: 1px dashed black; padding: 2px;">I</span>	Measurement angle (aperture stop position)	Integer; 1 digit	0: 1.0° 1: 0.2° 2: 0.1°
Explanation			
Reads the status of the instrument's aperture stop, which determines the measurement angle.			
Error-check codes			
Code	Meaning		
OK00	Normal completion		
ER00	Invalid command string/number of parameters		
ER83	Measurement angle abnormality		

UCCS (User Calibration Channel Select)			
Function			
Selects the user calibration channel to use.			
Input/Output Format			
PC		CS-2000/CS-2000A	
"UCCS,"1" + Delimiter code		⇒	
		⇐ "Error-check code" + Delimiter code	
Command Parameters			
	Meaning	Type	Details/range
1	User calibration channel	Integer; 2 digits	00 to 10 00: Konica Minolta calibration standard (no compensation) 1 to 10: User calibration channel to use
Response Parameters			
	Meaning	Type	Details/range
Explanation			
Sets the user calibration channel to use when determining measurement data. If set to a channel other than 00, user calibration data must already be stored in the instrument's internal memory for that channel.			
Error-check codes			
Code	Meaning		
OK00	Normal completion		
ER00	Invalid command string/number of parameters		
ER05	No compensation values (user calibration values).		
ER17	Parameter error Input user calibration channel is outside the setting range of 00 to 10.		
ER30	Instrument internal memory error		

UCCR ( <u>U</u> ser <u>C</u> alibration <u>C</u> hannel <u>R</u> ead)			
Function			
Reads the number of the currently selected user calibration channel.			
Input/Output Format			
PC		CS-2000/CS-2000A	
"UCCR" + <span style="border: 1px dashed black; padding: 2px;">Delimiter code</span>		⇒	
		⇐	" <span style="border: 1px dashed black; padding: 2px;">Error-check code</span> <span style="border: 1px dashed black; padding: 2px;">I</span> " + <span style="border: 1px dashed black; padding: 2px;">Delimiter</span> <span style="border: 1px dashed black; padding: 2px;">code</span>
Command Parameters			
	Meaning	Type	Details/range
Response Parameters			
	Meaning	Type	Details/range
<span style="border: 1px dashed black; padding: 2px;">I</span>	User calibration channel	Integer; 2 digits	00 to 10 00: Konica Minolta calibration standard (no compensation) 01 to 10: Currently selected user calibration channel • For single-digit channels, "0" will be added before the digit.
Explanation			
Reads the currently selected user calibration channel.			
Error-check codes			
Code	Meaning		
OK00	Normal completion		

UCPS (User Calibration Parameter Set)			
Function			
Sets user calibration parameters.			
Input/Output Format			
PC		CS-2000/CS-2000A	
Setting user calibration parameters is a three-step process:			
Step 1: Send user calibration coefficients. (Repeat for each wavelength.)			
"UCPS, [1], [2], [3], [4]" + <div>Delimiter code</div> ⇒			
Step 2: Send user calibration ID name.			
"UCPS, 2, [2], [5]" + <div>Delimiter code</div> ⇒			
Step 3: Finalize data and write to instrument internal memory.			
"UCPS, 3" + <div>Delimiter code</div> ⇒			
⇐ "Error-check code" + <div>Delimiter code</div>			
Command Parameters			
	Meaning	Type	Details/range
[1]	User calibration type	Integer; 1 digit	0: Wavelength correction 1: Level compensation
[2]	User calibration channel	Integer; 2 digits	1 to 10 User calibration channel to set data for.
[3]	Wavelength number	Integer; 3 digits	000 to 400 (380nm to 780nm) Wavelength for which calibration data will be written "000" = 380nm, "001" = 381nm, ... "400" = 780nm.
[4]	When Command Parameter [1] = "0" (Wavelength correction): Wavelength correction factor	Hex format	The corrected wavelength should be input. For example, to set 401nm as the corrected wavelength for 400nm, input "401" in hexadecimal format. Range: Nominal wavelength ±2nm
	When Command Parameter [1] = "1" (Level compensation): Level compensation factor	Hex format	For level compensation, the compensation factor should be input as absolute value, not percentage. (For example, 10% should be written as 0.1) Range: 0.001 to 1000
[5]	Calibration channel ID name	String; 10 char.	Alphanumeric string Length: 10 characters (if name is less than 10 characters, add spaces to achieve 10 characters) Refer to section 7: Characters.
Response Parameters			
	Meaning	Type	Details/range
Explanation			
Sets user calibration data and calibration channel ID. The procedure is a 3-step procedure; Data are not finalized and written to the instrument internal memory until step 3 is completed.			
Error-check codes			
Code	Meaning		
OK00	Normal completion		

ER00	Invalid command string/number of parameters
ER17	Parameter error Parameter set to value outside setting range
ER30	Instrument internal memory error

UCPR (User Calibration Parameter Read)			
Function			
Reads user calibration parameters for the specified user calibration channel.			
Input/Output Format			
PC		CS-2000/CS-2000A	
"UCPR,[1],[2] (,[3] )" + <span>Delimiter</span>		⇒	
<span>code</span>		⇐ " <span>Error-check code</span> ,[1]" + <span>Delimiter</span>	
		<span>code</span>	
Command Parameters			
	Meaning	Type	Details/range
[1]	User calibration data requested	Integer; 1 digit	0: Wavelength correction 1: Level compensation 2: Calibration channel ID name
[2]	User calibration channel	Integer; 2 digits	01 to 10 User calibration channel to read data from
[3]	When Command Parameter [1] = "0" (Wavelength correction) or "1" (Level compensation): Wavelength number	Integer; 3 digits	000 to 400 (380nm to 780nm) Wavelength for which calibration data will be read. "000" = 380nm, "001" = 381nm, ... "400" = 780nm.
	When Command Parameter [1] = "2" (Calibration channel ID name): Do not input		
Response Parameters			
	Meaning	Type	Details/range
[1]	When Command Parameter [1] = "0" (Wavelength correction): Wavelength correction factor	Hex format	The corrected wavelength is output.
	When Command Parameter [1] = "1" (Level compensation): Level compensation factor	Hex format	The compensation factor is output as absolute value, not percentage. (For example, 10% is output as 0.1) Range: 0.001 to 1000
	When Command Parameter [1] = "2" (Calibration channel ID name): Calibration channel ID name	String; 10 char.	Alphanumeric string Length: 10 characters (if name is less than 10 characters, additional spaces will be used to achieve 10 characters)
Explanation			
Reads user calibration data or ID name.			
Error-check codes			
Code	Meaning		
OK00	Normal completion		
ER00	Invalid command string/number of parameters		
ER05	No compensation values (user calibration values)		
ER17	Parameter error Parameter set to value outside setting range		
ER30	Instrument internal memory error		

UCCD (User Calibration Channel Delete)			
Function			
Deletes data stored in the specified user calibration channel.			
Input/Output Format			
PC		CS-2000/CS-2000A	
"UCCD," + <span>1</span> + <span>Delimiter code</span>		⇒	
		⇐ " <span>Error-check code</span> " + <span>Delimiter code</span>	
Command Parameters			
	Meaning	Type	Details/range
<span>1</span>	User calibration channel	Integer; 2 digits	1 to 10 User calibration channel to delete data from
Response Parameters			
	Meaning	Type	Details/range
Explanation			
Deletes data from the specified user calibration channel.			
Error-check codes			
Code	Meaning		
OK00	Normal completion		
ER00	Invalid command string/number of parameters		
ER17	Parameter error Calibration channel set to a value outside the range of 1 to 10.		
ER30	Instrument internal memory error		



LNSS (Close-up Lens Status Select)			
Function			
Selects the status of whether or not the optional close-up lens is attached.			
Input/Output Format			
PC		CS-2000/CS-2000A	
"LNSS," + <span style="border: 1px solid black; padding: 0 5px;">1</span> + <span style="border: 1px dashed black; padding: 0 5px;">Delimiter code</span>		⇒	
		⇐ " <span style="border: 1px dashed black; padding: 0 5px;">Error-check code</span> " + <span style="border: 1px dashed black; padding: 0 5px;">Delimiter code</span>	
Command Parameters			
	Meaning	Type	Details/range
<span style="border: 1px solid black; padding: 0 5px;">1</span>	Close-up lens status	Integer; 1 digit	0: None 1: Attached
Response Parameters			
	Meaning	Type	Details/range
Explanation			
<p>Sets whether or not the optional close-up lens is attached to the instrument.</p> <p>The CS-2000/CS-2000A does not automatically recognize whether or not a close-up lens is attached to the instrument. It is therefore necessary to use this command to tell the instrument whether or not a close-up lens is attached.</p> <p>When a close-up lens is attached, the lens compensation factors must be set in the instrument's memory in advance.</p>			
Error-check codes			
Code	Meaning		
OK00	Normal completion		
ER00	Invalid command string/number of parameters		
ER05	No compensation values (lens compensation factors)		
ER17	Parameter error.		
ER30	Instrument internal memory error		
ER83	Measurement angle abnormality		

LNSR (Close-up <u>L</u> ens <u>S</u> tatus <u>R</u> ead)			
Function			
Reads the status of whether or not the optional close-up lens is attached.			
Input/Output Format			
PC		CS-2000/CS-2000A	
"LNSR" + <span style="border: 1px dashed black; padding: 2px;">Delimiter code</span>		⇒	
		⇐	" <span style="border: 1px dashed black; padding: 2px;">Error-check code</span> <span style="border: 1px dashed black; padding: 2px;">I</span> " + <span style="border: 1px dashed black; padding: 2px;">Delimiter</span> <span style="border: 1px dashed black; padding: 2px;">code</span>
Command Parameters			
	Meaning	Type	Details/range
Response Parameters			
	Meaning	Type	Details/range
<span style="border: 1px dashed black; padding: 2px;">I</span>	Close-up lens status	Integer; 1 digit	0: None 1: Attached
Explanation			
Reads whether or not the optional close-up lens is attached to the instrument.			
Error-check codes			
Code	Meaning		
OK00	Normal completion		
ER00	Invalid command string/number of parameters		

ALFS (Attachment Lens Compensation Factor Set)

Function

Sets compensation factors for optional attachment lens (closeup lens).

Input/Output Format

PC

CS-2000/CS-2000A

Setting lens compensation factors is a two-step process:

Step 1: Send lens compensation factors. (Repeat for each wavelength.)

"ALFS, [1], [2], [3]" + [Delimiter code]⇒

Step 2: Finalize data and write to instrument internal memory.

"ALFS, 3" + [Delimiter code]⇒

⇐ "[Error-check code]" + [Delimiter code]

Command Parameters

	Meaning	Type	Details/range
[1]	Measurement angle	Integer; 1 digit	0: 1° 1: 0.2° 2: 0.1°
[2]	Wavelength number	Integer; 3 digits	000 to 400 (380nm to 780nm) Wavelength for which compensation factor will be written. "000" = 380nm, "001" = 381nm, ... "400" = 780nm.
[3]	Compensation factor	Hex format	The compensation factor should be input as absolute value, not percentage. (For example, 10% should be written as 0.1) Range: 0 to 1

Response Parameters

	Meaning	Type	Details/range

Explanation

Sets compensation factor when using optional close-up attachment lens. The factors which should be used are included with the optional close-up lens at the time of purchase.  
Data are not finalized and written to the instrument internal memory until "ALFS, 3" is sent.

Error-check codes

Code	Meaning
OK00	Normal completion
ER00	Invalid command string/number of parameters
ER17	Parameter error
ER30	Instrument internal memory error



NDFS (External <u>N</u> D <u>F</u> ilter <u>S</u> elect)			
Function			
Selects which optional external ND filter (if any) is attached.			
Input/Output Format			
PC		CS-2000/CS-2000A	
"NDFS," + <span style="border: 1px solid black; padding: 0 5px;">1</span> + <span style="border: 1px solid black; padding: 0 5px;">Delimiter code</span>		⇒	
		⇐	" <span style="border: 1px solid black; padding: 0 5px;">Error-check code</span> " + <span style="border: 1px solid black; padding: 0 5px;">Delimiter</span> <span style="border: 1px solid black; padding: 0 5px;">code</span>
Command Parameters			
	Meaning	Type	Details/range
<span style="border: 1px solid black; padding: 0 5px;">1</span>	External ND filter	Integer; 1 digit	0: None 1: ND Filter 1/10 attached 2: ND Filter 1/100 attached
Response Parameters			
	Meaning	Type	Details/range
Explanation			
<p>Sets which (if any) optional external ND filter is attached to the instrument.</p> <p>The CS-2000/CS-2000A does not automatically recognize whether or not an optional external ND filter is attached to the instrument. It is therefore necessary to use this command to tell the instrument whether or not an ND filter is attached, and, if attached, which one is attached.</p> <p>Two optional ND filters are available: ND Filter 1/10 and ND Filter 1/100. The proper setting must be made with this command. If the setting does not correspond to the attached ND filter, accurate measurements cannot be performed.</p> <p>When a close-up lens is attached, the lens compensation factors must be set in the instrument's memory in advance.</p>			
Error-check codes			
Code	Meaning		
OK00	Normal completion		
ER00	Invalid command string/number of parameters		
ER05	No compensation values (lens compensation factors)		
ER17	Parameter error Value set to other than 0, 1, or 2.		
ER83	Measurement angle abnormality		

NDFR (External <u>N</u> D <u>F</u> ilter <u>R</u> ead)			
Function			
Reads which optional external ND filter (if any) is attached.			
Input/Output Format			
PC		CS-2000/CS-2000A	
"NDFR" + <span style="border: 1px solid black; padding: 2px;">Delimiter code</span>		⇒	
		⇐	" <span style="border: 1px solid black; padding: 2px;">Error-check code</span> , <span style="border: 1px solid black; padding: 2px;">I</span> " + <span style="border: 1px solid black; padding: 2px;">Delimiter</span> <span style="border: 1px solid black; padding: 2px;">code</span>
Command Parameters			
	Meaning	Type	Details/range
Response Parameters			
	Meaning	Type	Details/range
<span style="border: 1px solid black; padding: 2px;">I</span>	External ND filter	Integer; 1 digit	0: None 1: ND Filter 1/10 attached 2: ND Filter 1/100 attached
Explanation			
Reads which (if any) optional external ND filter is attached to the instrument.			
Error-check codes			
Code	Meaning		
OK00	Normal completion		
ER00	Invalid command string/number of parameters		

NFCS (ND Filter Compensation Factor Set)			
Function			
Sets compensation factors for an optional external ND filter attached to instrument.			
Input/Output Format			
PC		CS-2000/CS-2000A	
Setting compensation factors is a two-step process:			
Step 1: Send compensation factors. (Repeat for each wavelength.)			
"NFCS, [1], [2], [3], [4]" + <span style="border: 1px dashed black; padding: 0 5px;">Delimiter code</span> ⇒			
Step 2: Finalize data and write to instrument internal memory.			
"NFCS, 3" + <span style="border: 1px dashed black; padding: 0 5px;">Delimiter code</span> ⇒			
⇐ " <span style="border: 1px dashed black; padding: 0 5px;">Error-check code</span> " + <span style="border: 1px dashed black; padding: 0 5px;">Delimiter code</span>			
Command Parameters			
	Meaning	Type	Details/range
[1]	Measurement angle	Integer; 1 digit	0: 1° 1: 0.2° 2: 0.1°
[2]	External ND filter	Integer; 1 digit	1: ND Filter 1/10 2: ND Filter 1/100
[3]	Wavelength number	Integer; 3 digits	000 to 400 (380nm to 780nm) Wavelength for which compensation factor will be written. "000" = 380nm, "001" = 381nm, ... "400" = 780nm.
[4]	Compensation factor	Hex format	The compensation factor should be input as absolute value, not percentage. (For example, 10% should be written as 0.1) Range: 0 to 1
Response Parameters			
	Meaning	Type	Details/range
Explanation			
Sets compensation factor when using optional external ND filter lens. The factors which should be used are included with the optional ND filter at the time of purchase. Data are not finalized and written to the instrument internal memory until "NFCS , 3" is sent.			
Error-check codes			
Code	Meaning		
OK00	Normal completion		
ER00	Invalid command string/number of parameters		
ER17	Parameter error Parameter set to value outside setting range		
ER30	Instrument internal memory error		
ER83	Measurement angle abnormality		

NFCR (ND Filter Compensation Factor Read)			
Function			
Reads compensation factors for external ND filter attached to instrument.			
Input/Output Format			
PC		CS-2000/CS-2000A	
"NFCR, [1], [2], [3]" + [Delimiter code]		⇒	
		⇐ "[Error-check code], [1]" + [Delimiter code]	
Command Parameters			
	Meaning	Type	Details/range
[1]	Measurement angle	Integer; 1 digit	0: 1° 1: 0.2° 2: 0.1°
[2]	External ND filter number	Integer; 1 digit	1: ND1 2: ND2
[3]	Wavelength number	Integer; 3 digits	000 to 400 (380nm to 780nm) Wavelength for which compensation factor will be written.
Response Parameters			
	Meaning	Type	Details/range
[3]	Compensation factor	Hex format	Compensation factor stored in the instrument for the specified measurement angle and wavelength
Explanation			
Reads compensation factors stored in instrument for when external ND filter is attached.			
Error-check codes			
Code	Meaning		
OK00	Normal completion		
ER00	Invalid command string/number of parameters		
ER05	No compensation values (ND filter compensation factors)		
ER17	Parameter error Parameter set to value outside setting range		



MEAS (Measure)			
Function			
Performs measurement or cancels measurement in progress.			
Input/Output Format			
PC		CS-2000/CS-2000A	
(To perform measurement)			
"MEAS, [1]" + [Delimiter code]		⇒	(Instrument performs pre-measurement. Time required: About 1 to 10s)
		⇐	"[Error-check code], [1]" + [Delimiter code] (Instrument starts actual measurement for time indicated by Response Parameter 1.)
		⇐	"[Error-check code]" + [Delimiter code] (Measurement completed.)
(To cancel measurement in progress)			
"MEAS, [1]" + [Delimiter code]		⇒	(Measurement is canceled.)
		⇐	"[Error-check code]" + [Delimiter code]
Command Parameters			
	Meaning	Type	Details/range
[1]	Command parameter	Integer; 1 digit	0: Cancel measurement 1: Start measurement
Response Parameters			
	Meaning	Type	Details/range
[1]	Measurement time	Integer; 3 digits	002 to 242: 3-character string indicating measurement time in seconds (as determined by pre-measurement) from time of response.
Explanation			
<p>To perform measurement: Measurement process starts when "MEAS,1" is input. A pre-measurement is taken to determine the required measurement time (and notification of this time is sent from the CS-2000/CS-2000A to the PC) and then the actual measurement begins automatically. When measurement has been completed, the instrument returns an error-check code ("OK00" if measurement was completed successfully).</p> <p>To cancel a measurement in progress, "MEAS,0" can be input after the pre-measurement has been completed. No commands will be accepted during pre-measurement. During actual measurement, commands other than "MEAS,0" will result in a response of "ER00".</p>			
Error-check codes			
Code	Meaning		
OK00	Normal completion		
ER00	Invalid command string/number of parameters		
ER10	Over measurement range		

ER17	Parameter error: "MEAS,0" sent when measurement is not in progress. "MEAS,1" sent while measurement is already in progress.
ER51 ER52	Temperature abnormality
ER71	Outside synchronization signal range
ER83	Measurement angle abnormality

<b>MEDR (<u>M</u>easurement <u>D</u>ata <u>R</u>ead)</b>			
<b>Function</b>			
Reads measurement data from instrument.			
<b>Input/Output Format</b>			
PC		CS-2000/CS-2000A	
(For normal measurement)			
<p>"MEDR,[1],[2],[3]" + [Delimiter code]      ⇒</p> <p>  ⇐ "[Error-check code],[1]</p> <p>   (,[2],...)" + [Delimiter code]</p>			
<b>Command Parameters</b>			
	Meaning	Type	Details/range
[1]	Data mode	Integer; 1 digit	0: Measurement conditions 1: Spectral data 2: Colorimetric data
[2]	Data format	Integer; 1 digit	0: Alphanumeric 1: Hexadecimal
[3]	Data block number to read	Integer; Up to 3 digits	<u>When Command Parameter [1] = "0" (Measurement conditions):</u> 1 (fixed)  <u>When Command Parameter [1] = "1" (Spectral data):</u> Wavelength group to read: 1: 100 pieces of data from 380 to 479nm 2: 100 pieces of data from 480 to 579nm 3: 100 pieces of data from 580 to 679nm 4: 101 pieces of data from 680 to 780nm  <u>When Command Parameter [1] = "2" (Colorimetric data):</u> 00: All colorimetric data 01: X,Y,Z 02: x, y, Lv 03: u',v',Lv 04: T, Δuv, Lv 05: λd, Pe,Lv 11: X <sub>10</sub> ,Y <sub>10</sub> ,Z <sub>10</sub> 12: x <sub>10</sub> ,y <sub>10</sub> ,Lv <sub>10</sub> 13: u' <sub>10</sub> ,v' <sub>10</sub> ,Lv <sub>10</sub> 14: T <sub>10</sub> ,Δuv <sub>10</sub> ,Lv <sub>10</sub> 15: λd <sub>10</sub> ,Pe <sub>10</sub> ,Lv <sub>10</sub> 100: Le 101: Lv

Response Parameters			
	Meaning	Type	Details/range
<u>When Command Parameter [1] =0 (Measurement conditions):</u>			
[1]	Speed mode:	Integer; 1 digit	0: Normal 1: Fast 2: Multi Integ Normal 3: Manual 4: Multi Integ Fast • For CS-2000 with firmware ver. 1.01.0000 or earlier, "2" means "Multi Integ" and "4" will not be output.
[2]	Sync mode:	Integer; 1 digit	0: No sync 1: Internal sync 2: External sync
[3]	Integration time	Integer; 9 digits	Integration time in $\mu\text{sec}$ • If number of digits is fewer than 9, "0" will be added before value.
[4]	Internal ND filter:	Integer; 1 digit	0: Off                      1: On
[5]	Optional close-up lens:	Integer; 1 digit	0: None                      1: Attached (Setting stored on instrument. Not automatically detected.)
[6]	Optional external ND filter:	Integer; 1 digit	0: None 1: ND Filter 1/10 2: ND Filter 1/100 (Setting stored on instrument. Not automatically detected.)
[7]	Measurement angle:	Integer; 1 digit	0: 1°                      1: 0.2° 2: 0.1°
[8]	Calibration channel:	Integer; 2 digits	00 to 10 00:                      Konica Minolta calibration standard (no compensation) 01 to 10: User calibration channel
<u>When Command Parameter [1] =1 (Spectral data):</u>			
[1]	Block of spectral irradiance data		(For type and format, please refer to Section 5: Numerical output formats)
<u>When Command Parameter [1] =2 (Colorimetric data):</u>			
[1]	Selected block of colorimetric data		(For type and format, please refer to Section 5: Numerical output formats)
Explanation			
The most recent measured data are read from the instrument memory.			
The timing for clearing measurement data from the instrument memory differs depending on whether or not the measurement button is enabled.			
When measurement button is not enabled:			
Measurement data are cleared when the next measurement is started.			
When measurement button is disabled:			
Measurement data are cleared when reading of all 4 blocks of spectral data has been completed or when reading of any 1 set of colorimetric data has been completed.			
Error-check codes			
Code	Meaning		
OK00	Normal completion		

ER00	Invalid command string/number of parameters
ER02	Measurement in process
ER10	Over measurement range
ER17	Parameter error
ER20	No measurement data
ER51 ER52	Temperature abnormality
ER71	Outside synchronization signal range
ER83	Measurement angle abnormality

MSWE ( <u>M</u> easuring <u>S</u> witch <u>E</u> nable)			
Function			
Enables/disables the measuring button in remote mode.			
Input/Output Format			
PC		CS-2000/CS-2000A	
"MSWE," + <i>Delimiter code</i>		⇒	
		⇐ "Error-check code" + <i>Delimiter code</i>	
Command Parameters			
	Meaning	Type	Details/range
"1"	Measuring button status	Integer; 1 digit	0: Disabled 1: Enabled
Response Parameters			
	Meaning	Type	Details/range
Explanation			
Enables/disables the measuring button.			
Error-check codes			
Code	Meaning		
OK00	Normal completion		
ER00	Invalid command string/number of parameters		
ER17	Parameter error		

STDS ( <u>S</u> tore <u>D</u> ata <u>S</u> et)			
Function			
Stores current measurement data to memory number.			
Input/Output Format			
PC		CS-2000/CS-2000A	
"STDS," + <span style="border: 1px dashed black; padding: 0 5px;">1</span> + <span style="border: 1px dashed black; padding: 0 5px;">Delimiter code</span>		⇒	
		⇐	" <span style="border: 1px dashed black; padding: 0 5px;">Error-check code</span> " + <span style="border: 1px dashed black; padding: 0 5px;">Delimiter</span> <span style="border: 1px dashed black; padding: 0 5px;">code</span>
Command Parameters			
	Meaning	Type	Details/range
<span style="border: 1px dashed black; padding: 0 5px;">1</span>	Memory number to store data in	Integer; 2 digits	0 to 99
Response Parameters			
	Meaning	Type	Details/range
Explanation			
Copies the most recent measured data to the specified memory number. If data already exists in that memory number, the existing data will be overwritten.			
Error-check codes			
Code	Meaning		
OK00	Normal completion		
ER00	Invalid command string/number of parameters		
ER17	Parameter error Memory number set to a value outside the range of 00 to 99.		
ER20	No data No measurement data available for copying to memory number.		
ER30	Instrument internal memory error		

STDR (Stored Measurement Data Read)

Function

Reads stored measurement data from instrument.

Input/Output Format

PC

CS-2000/CS-2000A

(For normal measurement)

"STDR,[1],[2],[3],[4]" + Delimiter code ⇒  
⇐ "[Error-check code],[1]  
(,[2],...)" + Delimiter code

Command Parameters

	Meaning	Type	Details/range
[1]	Memory number to read data from	Integer; 2 digits	00 to 99 <ul style="list-style-type: none"><li>For single-digit numbers, add a "0" before the digit.</li></ul>
[2]	Data mode	Integer; 1 digit	0: Measurement conditions 1: Spectral data 2: Colorimetric data
[3]	Data format	Integer; 1 digit	0: Alphanumeric 1: Hexadecimal (IEEE floating point format: 4-byte big-endian hexadecimal string)
[4]	Data block number to read	Integer; Up to 3 digits	<div>When Command Parameter [2] = "0" (Measurement conditions): Must be "1".</div> <div>When Command Parameter [2] = "1" (Spectral data): Wavelength group to read: 1: 100 pieces of data from 380 to 479nm 2: 100 pieces of data from 480 to 579nm 3: 100 pieces of data from 580 to 679nm 4: 101 pieces of data from 680 to 780nm</div> <div>When Command Parameter [2] = "2" (Colorimetric data): 00: All colorimetric data 01: X,Y,Z 02: x, y, Lv 03: u',v',Lv 04: T, Δuv, Lv 05: λd, Pe,Lv 11: X<sub>10</sub>,Y<sub>10</sub>,Z<sub>10</sub> 12: x<sub>10</sub>,y<sub>10</sub>,Lv<sub>10</sub> 13: u'<sub>10</sub>,v'<sub>10</sub>,Lv<sub>10</sub> 14: T<sub>10</sub>,Δuv<sub>10</sub>,Lv<sub>10</sub> 15: λd<sub>10</sub>,Pe<sub>10</sub>,Lv<sub>10</sub> 100: Le 101: Lv</div>



Response Parameters			
	Meaning	Type	Details/range
<u>When Command Parameter [1] =0 (Measurement conditions):</u>			
[1]	Speed mode:	Integer; 1 digit	0: Normal 1: Fast 2: Multi Integ Normal 3: Manual 4: Multi Integ Fast • For CS-2000 with firmware ver. 1.01.0000 or earlier, "2" means "Multi Integ" and "4" will not be output.
[2]	Sync mode:	Integer; 1 digit	0: No sync 1: Internal sync 2: External sync
[3]	Internal ND filter:	Integer; 1 digit	0: Off                      1: On
[4]	Optional close-up lens:	Integer; 1 digit	0: None                      1: Attached (Current setting stored on instrument. Not automatically detected.)
[5]	Optional external ND filter:	Integer; 1 digit	0: None 1: ND Filter 1/10 2: ND Filter 1/100 (Current setting stored on instrument.)
[6]	Measurement angle:	Integer; 1 digit	0: 1°                      1: 0.2° 2: 0.1°
[7]	Calibration channel:	Integer; 2 digits	00 to 10 00:                      Konica Minolta calibration standard (no compensation) 01 to 10: User calibration channel to use
<u>When Command Parameter [1] =1 (Spectral data):</u>			
[7]	Block of spectral irradiance data		(For type and format, please refer to Section 5: Numerical output formats)
<u>When Command Parameter [1] =2 (Colorimetric data):</u>			
[7]	Selected block of colorimetric data		(For type and format, please refer to Section 5: Numerical output formats)
Explanation			
Reads data stored in the specified memory channel from the instrument.			
Error-check codes			
Code	Meaning		
OK00	Normal completion		
ER00	Invalid command string/number of parameters		
ER17	Parameter error Input parameter is outside setting range.		
ER20	No data No data are stored in the specified target channel.		

STDD ( <u>S</u> tore <u>D</u> ata <u>D</u> ele)			
Function			
Deletes data stored in the specified memory number.			
Input/Output Format			
PC		CS-2000/CS-2000A	
"STDD," + <span>[1]</span> + <span>Delimiter code</span>		⇒	
		⇐ " <span>Error-check code</span> " + <span>Delimiter code</span>	
Command Parameters			
	Meaning	Type	Details/range
<span>[1]</span>	Memory number to delete data from	Integer; 2 digits	00 to 99
Response Parameters			
	Meaning	Type	Details/range
Explanation			
Deletes data from the specified memory number.			
Error-check codes			
Code	Meaning		
OK00	Normal completion		
ER00	Invalid command string/number of parameters		
ER17	Parameter error Memory number set to a value outside the range of 00 to 99.		

STAD ( <u>S</u> tore <u>A</u> ll <u>D</u> ata <u>D</u> ele <del>t</del> e)			
Function			
Deletes data stored in all memory numbers.			
Input/Output Format			
PC		CS-2000/CS-2000A	
"STAD" + <span style="border: 1px dashed black; padding: 2px;">Delimiter code</span>		⇒	
		⇐	" <span style="border: 1px dashed black; padding: 2px;">Error-check code</span> " + <span style="border: 1px dashed black; padding: 2px;">Delimiter</span> <span style="border: 1px dashed black; padding: 2px;">code</span>
Command Parameters			
	Meaning	Type	Details/range
Response Parameters			
	Meaning	Type	Details/range
Explanation			
Deletes data from all memory numbers.			
Error-check codes			
Code	Meaning		
OK00	Normal completion		
ER00	Invalid command string/number of parameters		

TGSL ( <u>T</u> arget <u>N</u> umber <u>S</u> elect)			
Function			
Selects the target number.			
Input/Output Format			
PC		CS-2000/CS-2000A	
"TGSL," + <span style="border: 1px dashed black; padding: 0 5px;">1</span> + <span style="border: 1px dashed black; padding: 0 5px;">Delimiter code</span>		⇒	
		⇐	" <span style="border: 1px dashed black; padding: 0 5px;">Error-check code</span> " + <span style="border: 1px dashed black; padding: 0 5px;">Delimiter</span> <span style="border: 1px dashed black; padding: 0 5px;">code</span>
Command Parameters			
	Meaning	Type	Details/range
<span style="border: 1px dashed black; padding: 0 5px;">1</span>	Target number	Integer; 2 digits	01 to 20: Target number • If no target data have been set, the target number should be set to 01.
Response Parameters			
	Meaning	Type	Details/range
Explanation			
Selects the target number to use for displaying color difference . When a target number other than "00" is specified, target color data must already be stored in that the memory channel for that number.			
Error-check codes			
Code	Meaning		
OK00	Normal completion		
ER00	Invalid command string/number of parameters		
ER17	Parameter error		

TGSR (Target Number Selection Read)			
Function			
Reads the currently selected target number.			
Input/Output Format			
PC		CS-2000/CS-2000A	
"TGSR" + <span style="border: 1px dashed black; padding: 2px;">Delimiter code</span>		⇒	
		⇐	" <span style="border: 1px dashed black; padding: 2px;">Error-check code</span> , <span style="border: 1px dashed black; padding: 2px;">1</span> " + <span style="border: 1px dashed black; padding: 2px;">Delimiter code</span>
Command Parameters			
	Meaning	Type	Details/range
Response Parameters			
	Meaning	Type	Details/range
<span style="border: 1px dashed black; padding: 2px;">1</span>	Target number	Integer; 2 digits	01 to 20 If no target data have been set, 01 will be returned as the target number. <ul style="list-style-type: none"> <li>If the target number is a single digit, a "0" will be added in front of the number.</li> </ul>
Explanation			
Reads the currently selected target number to use for displaying color difference.			
Error-check codes			
Code	Meaning		
OK00	Normal completion		
ER00	Invalid command string/number of parameters		

TGDS (Target Data Set)			
Function			
Sets target data in the specified target number channel of the instrument.			
Input/Output Format			
PC		CS-2000/CS-2000A	
(For writing target data)			
"TGDS, [1], [2], [3], [4], [5]" + [Delimiter]		⇒	[code]
		⇐	"[Error-check code]" + [Delimiter] [code]
(For writing target ID name: Data format=0, Data mode =2)			
"TGDS, 0, 2, [3], [6]" + [Delimiter code]		⇒	
		⇐	"[Error-check code]" + [Delimiter] [code]
Command Parameters			
	Meaning	Type	Details/range
[1]	Data format	Integer; 1 digit	0: Alphanumeric 1: Hexadecimal
[2]	Data mode	Integer; 1 digit	0: Spectral data 1: Colorimetric data 2: Target ID name
[3]	Target number	Integer; 2 digits	01 to 20
[4]	When Command Parameter [2] = "0" (Spectral data): Wavelength number	Integer; 3 digits	000 to 400 (380nm to 780nm) Wavelength for which calibration data will be read. "000" = 380nm, "001" = 381nm, ... "400" = 780nm.
	When Command Parameter [2] = "1" (Colorimetric data): Colorimetric data number	Integer; 2 digits	00: All colorimetric data 01: X, Y, Z 02: x, y, Lv 03: u', v', Lv 04: T, Δuv, Lv 05: λd, Pe, Lv 11: X10, Y10, Z10 12: x10, y10, Lv10 13: u'10, v'10, Lv10 14: T10, Δuv10, Lv10 15: λd10, Pe10, Lv10 100: Le 101: Lv

5	When Command Parameter [2] = "0" (Spectral data): Spectral irradiance		(For type and format, please refer to Section 6: Numerical input formats)
	When Command Parameter [2] = "1" (Colorimetric data): Colorimetric data number		(For type and format, please refer to Section 6: Numerical input formats)
6	Target ID name	String; 10 char.	Alphanumeric string Length: 10 characters (if name is less than 10 characters, additional spaces will be used to achieve 10 characters) (Refer to Section 7: Acceptable characters)
<b>Response Parameters</b>			
	Meaning	Type	Details/range
<b>Explanation</b>			
<p>Sets target spectral irradiance, target colorimetric data, or target ID name in the selected target number channel of the instrument.</p> <p>For target spectral irradiance, data for all wavelengths must be written. If only part of the spectral irradiance data are input, the data are not stored in instrument internal memory.</p> <p>When spectral irradiance data are set, the colorimetric data are calculated and stored.</p>			
<b>Error-check codes</b>			
Code	Meaning		
OK00	Normal completion		
ER00	Invalid command string/number of parameters		
ER17	Parameter error This error is also output if an error occurs during calculation of colorimetric values.		
ER30	Instrument internal memory error		

TGDR (Target Data Read)			
Function			
Reads target data stored in the specified target number channel of the instrument.			
Input/Output Format			
PC		CS-2000/CS-2000A	
(For reading target data)			
"TGDS, [1],[2],[3],[4]" + [Delimiter code]		⇒	
		⇐ "[Error-check code],[1] (,[2],...)" + [Delimiter code]	
(For reading target ID name: Data format=0, Data mode =2)			
"TGDS, [1],[2]" + [Delimiter code]		⇒	
		⇐ "[Error-check code],[1]" + [Delimiter code]	
Command Parameters			
	Meaning	Type	Details/range
[1]	Target number	Integer; 2 digits	1 to 20
[2]	Data mode	Integer; 1 digit	0: Measurement conditions 1: Spectral data 2: Colorimetric data 3: Target ID name
[3]	Data format	Integer; 1 digit	0: Alphanumeric 1: Hexadecimal
[4]	Data block number to read	Integer; 3 digits	<div>When Command Parameter [2] = "0" (Measurement conditions): 1 (fixed)</div> <div>When Command Parameter [2] = "1" (Spectral data): Wavelength group to read: 01: 100 pieces of data from 380 to 479nm 02: 100 pieces of data from 480 to 579nm 03: 100 pieces of data from 580 to 679nm 04: 101 pieces of data from 680 to 780nm</div> <div>When Command Parameter [2] = "2" (Colorimetric data): Colorimetric data number 00: All colorimetric data 01: X,Y,Z 02: x, y, Lv 03: u',v',Lv 04: T, Δuv, Lv 05: λd, Pe,Lv 11: X10,Y10,Z10 12: x10, y10, Lv10 13: u'10,v'10,Lv10 14: T10, Δuv10, Lv10 15: λd10, Pe10,Lv10 100: Le</div>



			101: Lv
<b>Response Parameters</b>			
	Meaning	Type	Details/range
<u>When Command Parameter [2] =0 (Measurement conditions):</u>			
[1]	Speed mode:	Integer; 1 digit	0: Normal 1: Fast 2: Multi Integ Normal 3: Manual 4: Multi Integ Fast • For CS-2000 with firmware ver. 1.01.0000 or earlier, "2" means "Multi Integ" and "4" will not be output.
[2]	Sync mode:	Integer; 1 digit	0: No sync 1: Internal sync 2: External sync
[3]	Integration time	Integer; 9 digits	Integration time in $\mu\text{sec}$ • If number of digits is fewer than 9, "0" will be added before value.
[4]	Internal ND filter:	Integer; 1 digit	0: Off                      1: On
[5]	Optional close-up lens:	Integer; 1 digit	0: None                      1: Attached (Setting stored on instrument. Not automatically detected.)
[6]	Optional external ND filter:	Integer; 1 digit	0: None 1: ND Filter 1/10 2: ND Filter 1/100 (Setting stored on instrument. Not automatically detected.)
[7]	Measurement angle:	Integer; 1 digit	0: 1°                      1: 0.2° 2: 0.1°
[9]	Calibration channel:	Integer; 2 digits	00 to 10 00: Konica Minolta calibration standard (no compensation) 1 to 10: User calibration channel
<u>When Command Parameter [2] =1 (Spectral data):</u>			
[1]	Block of spectral irradiance data		(For type and format, please refer to Section 5: Numerical output formats)
<u>When Command Parameter [2] =2 (Colorimetric data):</u>			
[1]	Selected block of colorimetric data		(For type and format, please refer to Section 5: Numerical output formats)
<u>When Command Parameter [2] =3 (Target ID name):</u>			
[1]	Target ID name	String; 10 char.	Alphanumeric string Length: 10 characters (if name is less than 10 characters, additional spaces will be used to achieve 10 characters)
<b>Explanation</b>			

Reads target data stored in the specified target number channel of the instrument.  
Data must be already stored in the specified target channel

**Error-check codes**

Code	Meaning
OK00	Normal completion
ER00	Invalid command string/number of parameters
ER17	Parameter error A parameter was set outside the setting range.
ER20	No data are stored in the specified memory channel.

TGDD (Target Data Delete)			
Function			
Deletes data for the specified target number channel.			
Input/Output Format			
PC		CS-2000/CS-2000A	
"TGDD," + <i>Delimiter code</i>		⇒	
		⇐ " <i>Error-check code</i> " + <i>Delimiter code</i>	
Command Parameters			
	Meaning	Type	Details/range
<i>I</i>	Target number to delete	Integer; 2 digits	1 to 20
Response Parameters			
	Meaning	Type	Details/range
Explanation			
Deletes target data and target ID name from the specified target number.			
Error-check codes			
Code	Meaning		
OK00	Normal completion		
ER00	Invalid command string/number of parameters		
ER17	Parameter error		
ER30	Instrument internal memory error		

TGAD (Target All Data Delete)			
Function			
Deletes data stored in all target number channels.			
Input/Output Format			
PC		CS-2000/CS-2000A	
"TGAD" + <span style="border: 1px dashed black; padding: 2px;">Delimiter code</span>		⇒	
		⇐	" <span style="border: 1px dashed black; padding: 2px;">Error-check code</span> " + <span style="border: 1px dashed black; padding: 2px;">Delimiter</span> <span style="border: 1px dashed black; padding: 2px;">code</span>
Command Parameters			
	Meaning	Type	Details/range
Response Parameters			
	Meaning	Type	Details/range
Explanation			
Deletes target data and target ID name from all target numbers.			
Error-check codes			
Code	Meaning		
OK00	Normal completion		
ER00	Invalid command string/number of parameters		

BALS (Backlight Control Set)			
Function			
Sets control of external display backlight during measurement.			
Input/Output Format			
PC		CS-2000/CS-2000A	
"BALS, [1], [2]" + <span>Delimiter code</span>		⇒	
		⇐ " <span>Error-check code</span> " + <span>Delimiter code</span>	
Command Parameters			
	Meaning	Type	Details/range
[1]	External display backlight normally on/off	Integer; 1 digit	0: Off 1: On
[2]	External display backlight on/off during measurement	Integer; 1 digit	0: Off 1: On
Response Parameters			
	Meaning	Type	Details/range
Explanation			
Controls whether the backlight for the external display (LCD) on the rear of the instrument is on or off normally and during measurement. The results of the combined settings is shown below:			
		Command Parameter 2	
		0: Off	1: On
Command Parameter 1	0: Off	Measurement values shown:	Off
		During measurement:	Off
	1: On	Measurement values shown:	On
		During measurement:	On
Error-check codes			
Code	Meaning		
OK00	Normal completion		
ER00	Invalid command string/number of parameters		
ER17	Parameter error (Parameter set to a value outside the setting range of 0 or 1.)		

<b>BALR (Backlight Control Read)</b>			
<b>Function</b>			
Reads the setting for control of external display backlight during measurement.			
<b>Input/Output Format</b>			
PC		CS-2000/CS-2000A	
"BALR" + <span style="border: 1px dashed black; padding: 2px;">Delimiter code</span>		⇒	
		⇐	" <span style="border: 1px dashed black; padding: 2px;">Error-check code</span> , <span style="border: 1px dashed black; padding: 2px;">1</span> , <span style="border: 1px dashed black; padding: 2px;">2</span> "
		+	<span style="border: 1px dashed black; padding: 2px;">Delimiter code</span>
<b>Command Parameters</b>			
	Meaning	Type	Details/range
<b>Response Parameters</b>			
	Meaning	Type	Details/range
<span style="border: 1px dashed black; padding: 2px;">1</span>	External display backlight normally on/off	Integer; 1 digit	0: Off 1: On
<span style="border: 1px dashed black; padding: 2px;">2</span>	External display backlight on/off during measurement	Integer; 1 digit	0: Off 1: On
<b>Explanation</b>			
Reads the setting for control of whether the backlight for the external display (LCD) on the rear of the instrument is on or off normally and during measurement. For further information, please see the <u>BALS</u> command.			
<b>Error-check codes</b>			
Code	Meaning		
OK00	Normal completion		
ER00	Invalid command string/number of parameters		

CSMS ( <u>C</u> olor <u>S</u> pace <u>M</u> ode <u>S</u> et)			
Function			
Sets the color space to be displayed on the instrument's LCD.			
Input/Output Format			
PC		CS-2000/CS-2000A	
"CSMS," + <span style="border: 1px dashed black; padding: 0 5px;">1</span> + <span style="border: 1px dashed black; padding: 0 5px;">Delimiter code</span>		⇒	
		⇐	" <span style="border: 1px dashed black; padding: 0 5px;">Error-check code</span> " + <span style="border: 1px dashed black; padding: 0 5px;">Delimiter</span> <span style="border: 1px dashed black; padding: 0 5px;">code</span>
Command Parameters			
	Meaning	Type	Details/range
<span style="border: 1px dashed black; padding: 0 5px;">1</span>	Color space mode	Integer; 1 digit	0: L <sub>v</sub> , x, y 1: L <sub>v</sub> , u', v' 2: L <sub>v</sub> , T, Δuv 3: X, Y, Z 4: λ <sub>d</sub> , P <sub>e</sub> 5: Spectral graph
Response Parameters			
	Meaning	Type	Details/range
Explanation			
Sets the color space to be displayed on the LCD at the rear of the instrument.			
Error-check codes			
Code	Meaning		
OK00	Normal completion		
ER00	Invalid command string/number of parameters		
ER17	Parameter error Color space mode set to a value outside of the range 0 to 5 (inclusive)		

CSMR (Color Space Mode Read)			
Function			
Reads the color space to be displayed on the instrument's LCD.			
Input/Output Format			
PC		CS-2000/CS-2000A	
"CSMR" + <span style="border: 1px dashed black; padding: 2px;">Delimiter code</span>		⇒	
		⇐	" <span style="border: 1px dashed black; padding: 2px;">Error-check code</span> , <span style="border: 1px dashed black; padding: 2px;">1</span> " + <span style="border: 1px dashed black; padding: 2px;">Delimiter code</span>
Command Parameters			
	Meaning	Type	Details/range
Response Parameters			
	Meaning	Type	Details/range
<span style="border: 1px dashed black; padding: 2px;">1</span>	Color space mode	Integer; 1 digit	0: L <sub>v</sub> , x, y 1: L <sub>v</sub> , u', v' 2: L <sub>v</sub> , T, Δuv 3: X, Y, Z 4: λd, Pe 5: Spectral graph
Explanation			
Reads the color space to be displayed on the LCD at the rear of the instrument.			
Error-check codes			
Code	Meaning		
OK00	Normal completion		
ER00	Invalid command string/number of parameters		



DIMS (Display Mode Set)			
Function			
Sets the display mode (absolute or difference) for the instrument's LCD.			
Input/Output Format			
PC		CS-2000/CS-2000A	
"DIMS," + <span>1</span> + <span>Delimiter code</span>		⇒	
		⇐ " <span>Error-check code</span> " + <span>Delimiter code</span>	
Command Parameters			
	Meaning	Type	Details/range
<span>1</span>	Display mode	Integer; 1 digit	0: Absolute data display 1: Difference data display
Response Parameters			
	Meaning	Type	Details/range
Explanation			
Sets the display mode to show absolute colorimetric values (ABS) or colorimetric difference values (DIFF) in the LCD at the rear of the instrument.			
Error-check codes			
Code	Meaning		
OK00	Normal completion		
ER00	Invalid command string/number of parameters		
ER17	Parameter error		

DIMR ( <u>D</u> isplay <u>M</u> ode <u>R</u> ead)			
Function			
Reads the display mode (absolute or difference) for the instrument's LCD.			
Input/Output Format			
PC		CS-2000/CS-2000A	
"DIMR" + <span style="border: 1px dashed black; padding: 0 5px;">Delimiter code</span>		⇒	
		⇐	" <span style="border: 1px dashed black; padding: 0 5px;">Error-check code</span> , <span style="border: 1px dashed black; padding: 0 2px;">1</span> " + <span style="border: 1px dashed black; padding: 0 5px;">Delimiter</span> <span style="border: 1px dashed black; padding: 0 5px;">code</span>
Command Parameters			
	Meaning	Type	Details/range
Response Parameters			
	Meaning	Type	Details/range
<span style="border: 1px dashed black; padding: 0 2px;">1</span>	Display mode	Integer; 1 digit	0: Absolute data display 1: Difference data display
Explanation			
Reads the display mode indicating whether absolute colorimetric values (ABS) or colorimetric difference values (DIFF) are shown in the LCD at the rear of the instrument.			
Error-check codes			
Code	Meaning		
OK00	Normal completion		
ER00	Invalid command string/number of parameters		

OBSS ( <u>O</u> bserver <u>S</u> et)			
Function			
Sets the observer mode.			
Input/Output Format			
PC		CS-2000/CS-2000A	
"OBSS," + <span>1</span> + <span>Delimiter code</span>		⇒	
		⇐ " <span>Error-check code</span> " + <span>Delimiter code</span>	
Command Parameters			
	Meaning	Type	Details/range
<span>1</span>	Observer	Integer; 1 digit	0: 2° Standard Observer (2° OBS) 1: 10° Standard Observer (10° OBS)
Response Parameters			
	Meaning	Type	Details/range
Explanation			
Sets the CIE observer mode to 2° Standard Observer (2° OBS) or 10° Standard Observer (10° OBS) for colorimetric calculations.			
Error-check codes			
Code	Meaning		
OK00	Normal completion		
ER00	Invalid command string/number of parameters		
ER17	Parameter error		

OBSR ( <u>O</u> bserver <u>R</u> ead)			
Function			
Reads the observer mode.			
Input/Output Format			
PC		CS-2000/CS-2000A	
"OBSR" + <span style="border: 1px dashed black; padding: 2px;">Delimiter code</span>		⇒	
		⇐	" <span style="border: 1px dashed black; padding: 2px;">Error-check code</span> , <span style="border: 1px dashed black; padding: 2px;">i</span> " + <span style="border: 1px dashed black; padding: 2px;">Delimiter code</span>
Command Parameters			
	Meaning	Type	Details/range
Response Parameters			
	Meaning	Type	Details/range
<span style="border: 1px dashed black; padding: 2px;">i</span>	Observer		0: 2° Standard Observer (2° OBS) 1: 10° Standard Observer (10° OBS)
Explanation			
Reads the CIE observer mode for colorimetric calculations: 2° Standard Observer (2° OBS) or 10° Standard Observer (10° OBS).			
Error-check codes			
Code	Meaning		
OK00	Normal completion		
ER00	Invalid command string/number of parameters		

DTCR ( <u>D</u> ate/ <u>T</u> ime of <u>C</u> alibration <u>R</u> ead)		
Function		
Reads the date and time of factory calibration.		
Input/Output Format		
PC		CS-2000/CS-2000A
"DTCR" + <span style="border: 1px dashed black; padding: 2px;">Delimiter code</span> ⇨ ⇐ " <span style="border: 1px dashed black; padding: 2px;">Error-check code</span> , <span style="border: 1px dashed black; padding: 2px;">1</span> , <span style="border: 1px dashed black; padding: 2px;">2</span> " + <span style="border: 1px dashed black; padding: 2px;">Delimiter code</span>		
Command Parameters		
	Meaning	Details/range
Response Parameters		
	Meaning	Details/range
<span style="border: 1px dashed black; padding: 2px;">1</span>	Date of factory calibration	8 characters indicating year (4 characters), month (2 characters), and day (2 characters) For example, "20070201" means February 1, 2007.
<span style="border: 1px dashed black; padding: 2px;">2</span>	Time of factory calibration	6 characters indicating hour (2 characters; 24-hour clock), minute (2 characters), and second (2 characters). For example, "235607" indicates 23:56:07 (11:56:07 PM)
Explanation		
Reads the date and time of factory calibration.		
Error-check codes		
Code	Meaning	
OK00	Normal completion	
ER00	Invalid command string/number of parameters	

## 4 Error-check codes

Code	Meaning
OK00	Normal completion <ul style="list-style-type: none"> <li>Received command was processed normally</li> </ul>
ER00	<ul style="list-style-type: none"> <li>Received command string was not a valid command.</li> <li>Number of parameters received were incorrect for the command.</li> </ul>
ER02	Measurement in process <ul style="list-style-type: none"> <li>Received command cannot be processed because instrument is currently taking a measurement.</li> </ul>
ER05	No compensation values <ul style="list-style-type: none"> <li>There are no user calibration values in the specified calibration channel.</li> <li>There are no attachment lens compensation values in memory for the specified measurement angle.</li> <li>There are no ND filter compensation values in memory for the specified ND filter.</li> </ul>
ER10	Over measurement range <ul style="list-style-type: none"> <li>The luminance of the measurement subject exceeds the instrument's luminance measuring range.</li> <li>The flicker of the measurement subject is too large.</li> </ul>
ER17	Parameter error <ul style="list-style-type: none"> <li>The input parameter is outside the specified input range for the parameter (numerical value range or number of characters).</li> </ul>
ER20	No data <ul style="list-style-type: none"> <li>There are no measurement data in the instrument's memory buffer.</li> <li>There are no measurement data in the instrument's specified memory channel.</li> <li>There are no target data in the instruments specified target number channel.</li> </ul>
ER30	Instrument internal memory error <ul style="list-style-type: none"> <li>An error occurred while reading from or writing to the instrument's internal memory.</li> </ul>
ER51 ER52	Temperature abnormality <ul style="list-style-type: none"> <li>Ambient temperature during measurement is too high, causing the internal temperature of the instrument's sensor to become abnormal.</li> </ul>
ER71	Outside synchronization signal range <ul style="list-style-type: none"> <li>When instrument is set for external sync, vertical synchronization signal could not be detected.</li> <li>When instrument is set for external sync, vertical synchronization signal was below 20Hz or over 200Hz.</li> </ul>
ER81	Shutter operation abnormality <ul style="list-style-type: none"> <li>An abnormality occurred in the operation of the instrument's internal shutter mechanism.</li> </ul>
ER82	Internal ND filter operation malfunction <ul style="list-style-type: none"> <li>An abnormality occurred in the operation of the instrument's internal ND filter mechanism.</li> </ul>
ER83	Measurement angle abnormality <ul style="list-style-type: none"> <li>Measurement was performed with the measuring angle selector not set to a normal position.</li> <li>Measuring angle selector was moved during measurement.</li> </ul>
ER84	Cooling fan abnormality <ul style="list-style-type: none"> <li>Cooling fan is stopped.</li> <li>An abnormality occurred in the cooling mechanism.</li> </ul>
ER99	Program abnormality <ul style="list-style-type: none"> <li>An abnormality other than those covered by other error-check codes has occurred.</li> </ul>

## 5 Numerical output formats

### 5.1 Spectral irradiance values

When the CS-2000/CS-2000A outputs spectral irradiance (spectral data), the data output consists of comma-delimited spectral irradiance data at 1nm intervals for the data block number specified by the "Data block number" parameter sent with the command.

The wavelength range and number of data for each data block number are as follows:

Data block number	Wavelength range	Number of data
1	380 to 479nm	100
2	480 to 579nm	100
3	580 to 679nm	100
4	680 to 780nm	101

Each spectral irradiance data is output in either alphanumeric or hexadecimal format, according to the setting of the "Data format" parameter sent with the command.

#### 5.1.1 Alphanumeric data

Alphanumeric data are output in the following format:

- When an error occurs during internal calculations by the instrument, the calculation error number will be output.

Status	Format	Comment
Normal	Exponential format: #.####.e±# Decimal places: 4 digits Exponent: 1 digit	
Calculation error	"-9.9999e9"	For CS-2000 firmware ver. 1.01.0000 or earlier, "0" will be output.

#### 5.1.2 Hexadecimal data

Hexadecimal data are output in the following format:

- When an error occurs during internal calculations by the instrument, the calculation error number will be output.

Status	Format	Comment
Normal	Hexadecimal format: #### IEEE floating point format: 4-byte big-endian hexadecimal string (4 characters)	
Calculation error	-9.9999e10 in hexadecimal format	For CS-2000 firmware ver. 1.01.0000 or earlier, "0" will be output.

## 5.2 Colorimetric values

When the CS-2000/CS-2000A outputs colorimetric data, the data output consists of comma-delimited colorimetric data values for the data block number specified by the "Data block number" parameter sent with the command.

The colorimetric values and number and order of output data for each data block number are as follows:

Data block number	Colorimetric values	Data output
0	All colorimetric data	" Le , Lv , X , Y , Z , x , y , u' , v' , T , Δuv , λd , Pe , X <sub>10</sub> , Y <sub>10</sub> , Z <sub>10</sub> , x <sub>10</sub> , y <sub>10</sub> , u' <sub>10</sub> , v' <sub>10</sub> , T <sub>10</sub> , Δuv <sub>10</sub> , λd <sub>10</sub> , Pe <sub>10</sub> "
1	X,Y,Z	" X , Y , Z "
2	x, y, Lv	" x , y , Lv "
3	u',v',Lv	" u' , v' , Lv "
4	T, Δuv, Lv	" T , Δuv , Lv "
5	λd, Pe,Lv	" λd , Pe , Lv "
11	X <sub>10</sub> ,Y <sub>10</sub> ,Z <sub>10</sub>	" X <sub>10</sub> , Y <sub>10</sub> , Z <sub>10</sub> "
12	x <sub>10</sub> , y <sub>10</sub> , Lv	" x <sub>10</sub> , y <sub>10</sub> , Lv "
13	u' <sub>10</sub> ,v' <sub>10</sub> ,Lv	" u' <sub>10</sub> , v' <sub>10</sub> , Lv "
14	T <sub>10</sub> , Δuv <sub>10</sub> , Lv	" T <sub>10</sub> , Δuv <sub>10</sub> , Lv "
15	λd <sub>10</sub> , Pe <sub>10</sub> ,Lv	" λd <sub>10</sub> , Pe <sub>10</sub> , Lv "
100	Le	" Le "
101	Lv	" Lv "

Each colorimetric value is output in either alphanumeric or hexadecimal format, according to the setting of the "Data format" parameter sent with the command.



## 5.2.1 Alphanumeric data

Alphanumeric data are output in the following format:

- When an error occurs during internal calculations by the instrument, the calculation error number will be output.

Value	Status	Format	Comment
(9 characters)			
Le	Normal	Exponential format: “#.####e±#” Decimal places: 4 digits Exponent: 1 digit	
	Le<0.0001e-9	“0.0000e-9”	
	Calculation error	“-9.9999e9”	For CS-2000 firmware ver. 1.01.0000 or earlier, “0” will be output.
(6 characters)			
Lv	0.0001<Lv<10 0,000	Fixed character count of 6 characters; Floating-point decimal: “#####” or “####.#” or ... or “#.#####”	
	Lv>100,000	Unsigned exponential format: “#.##e#” Decimal places: 2 digits Exponent: 1 digit	For CS-2000 firmware ver. 1.01.0000 or earlier, there is no output in unsigned exponential format
	Lv<0.0001	“0.0000”	
	Calculation error	“-9.9e9”	For CS-2000 firmware ver. 1.01.0000 or earlier, “0” will be output.
(9 characters)			
X Y Z X <sub>10</sub> Y <sub>10</sub> Z <sub>10</sub>	Normal	Exponential format: “#.####e±#” Decimal places: 4 digits Exponent: 1 digit	
	Value<0.0001e-9	“0.0000e-9”	
	Calculation error	“-9.9999e9”	For CS-2000 firmware ver. 1.01.0000 or earlier, “0” will be output.
(6 characters)			
x y x <sub>10</sub> y <sub>10</sub>	Normal	Unsigned fixed-point decimal Decimal places: 4 digits “0.####”	
	Calculation error	“-9.999”	For CS-2000 firmware ver. 1.01.0000 or earlier, “0” will be output.
(6 characters)			
u' v' u <sub>10</sub> ' v <sub>10</sub> '	Normal	Unsigned fixed-point decimal Decimal places: 4 digits “0.####”	Upper limit: 0.9999 Lower limit: 0.0000
	Calculation error	“-9.999”	For CS-2000 firmware ver. 1.01.0000 or earlier, “0” will be output.
(5 characters)			
T T <sub>10</sub>	Normal	Unsigned integer; 5 digits maximum “#####”	Upper limit: 99999 Lower limit: 0
	Calculation error	“-9999”	For CS-2000 firmware ver. 1.01.0000 or earlier, “0” will be output.

Value	Status	Format	Comment
$\Delta u v$ $\Delta u v_{10}$	(7 characters)		
	Normal	Signed fixed-point decimal Decimal places: 4 digits “ $\pm 0 . \#\#\#$ ”	Upper limit: 0.9999 Lower limit: -0.9999
	Calculation error	“ $-9 . 9999$ ”	For CS-2000 firmware ver. 1.01.0000 or earlier, “0” will be output.
$\lambda d$ Pe	(6 characters)		
	Less than 100,000	Fixed character count of 6 characters; Floating-point decimal: “#####” or “####.#” or .... or “#.####”	Upper limit: 9.999e9 Lower limit: 0.0000
	Calculation error	“ $-9 . 9e9$ ”	For CS-2000 firmware ver. 1.01.0000 or earlier, “0” will be output.

### 5.2.2 Hexadecimal data

Hexadecimal data are output in the following format:

- When an error occurs during internal calculations by the instrument, the calculation error number will be output.

Status	Format	Comment
Normal	Hexadecimal format: ##### IEEE floating point format: 4-byte big-endian hexadecimal string (4 characters)	
Calculation error	-9.9999e10 in hexadecimal format	For CS-2000 firmware ver. 1.01.0000 or earlier, "0" will be output.

## 6 Numerical input formats

### 6.1 Spectral irradiance values

The command “TGDS” can be used to input target spectral irradiance values to the instrument in either alphanumeric or hexadecimal format, according to the setting of the “Data format” parameter sent with the command.

- All 401 spectral irradiance values (for wavelengths from 380 to 780nm) must be input for the instrument to store the target in memory. If fewer than 401 values are input, the target data will not be stored.
- When spectral irradiance values are input, the colorimetric values that will be the target color are calculated by the instrument internally.

#### 6.1.1 Alphanumeric data

The format for inputting alphanumeric data for target values is not specified.

- The total length of the characters sent to the instrument (including command, commas, and delimiter) must not exceed 800 bytes.

#### 6.1.2 Hexadecimal data

The format for inputting hexadecimal data for target values is as follows:

Format
Hexadecimal format: #####
IEEE floating point format: 4-byte big-endian hexadecimal string (4 characters)

## 6.2 Colorimetric values

The command "TGDS" can be used to input target comma-delimited colorimetric values for the data block number specified by the "Data block number" parameter sent with the command.

The colorimetric values and number and order of required input data for each data block number are as follows:

Data block number	Colorimetric values	Data to be input
0	All colorimetric data	" Le , Lv , X , Y , Z , x , y , u' , v' , T , Δuv , λd , Pe , X <sub>10</sub> , Y <sub>10</sub> , Z <sub>10</sub> , x <sub>10</sub> , y <sub>10</sub> , u' <sub>10</sub> , v' <sub>10</sub> , T <sub>10</sub> , Δuv <sub>10</sub> , λd <sub>10</sub> , Pe <sub>10</sub> "
1	X,Y,Z	" X , Y , Z "
2	x, y, Lv	" x , y , Lv "
3	u',v',Lv	" u' , v' , Lv "
4	T, Δuv, Lv	" T , Δuv , Lv "
5	λd, Pe,Lv	" λd , Pe , Lv "
11	X <sub>10</sub> ,Y <sub>10</sub> ,Z <sub>10</sub>	" X <sub>10</sub> , Y <sub>10</sub> , Z <sub>10</sub> "
12	x <sub>10</sub> , y <sub>10</sub> , Lv	" x <sub>10</sub> , y <sub>10</sub> , Lv "
13	u' <sub>10</sub> ,v' <sub>10</sub> ,Lv	" u' <sub>10</sub> , v' <sub>10</sub> , Lv "
14	T <sub>10</sub> , Δuv <sub>10</sub> , Lv	" T <sub>10</sub> , Δuv <sub>10</sub> , Lv "
15	λd <sub>10</sub> , Pe <sub>10</sub> ,Lv	" λd <sub>10</sub> , Pe <sub>10</sub> , Lv "
100	Le	" Le "
101	Lv	" Lv "

- When either [X, Y, Z], [x, y, Lv], or [u', v', Lv] are input, of the colorimetric values that will be the target color, [X, Y, Z], [x, y], [u', v'], [T, Δuv] and [λd, Pe] will be recalculated by the instrument internally.
- When either [X<sub>10</sub>, Y<sub>10</sub>, Z<sub>10</sub>], [x<sub>10</sub>, y<sub>10</sub>, Lv], or [u'<sub>10</sub>, v'<sub>10</sub>, Lv] are input, of the colorimetric values that will be the target color, [X<sub>10</sub>, Y<sub>10</sub>, Z<sub>10</sub>], [x<sub>10</sub>, y<sub>10</sub>], [u'<sub>10</sub>, v'<sub>10</sub>], [T<sub>10</sub>, Δuv<sub>10</sub>] and [λd<sub>10</sub>, Pe<sub>10</sub>] will be recalculated by the instrument internally.

Each colorimetric value is output in either alphanumeric or hexadecimal format, according to the setting of the "Data format" parameter sent with the command.

### 6.2.1 Alphanumeric data

The format for inputting alphanumeric data for target values is not specified.

- The total length of the characters sent to the instrument (including command, commas, and delimiter) must not exceed 800 bytes.

### 6.2.2 Hexadecimal data

The format for inputting hexadecimal data for target values is as follows:

Format
Hexadecimal format: ####
IEEE floating point format: 4-byte big-endian hexadecimal string (4 characters)

## 7 Characters for data communication

The following characters can be used for data communication with the CS-2000/CS-2000A:

- Numbers “0” to “9”
- Upper-case alphabetic letters “A” to “Z”
- Lower-case alphabetic letters “a” to “z”
- Space character

## 8 Measurement flow

### 8.1 Communication flow for measurement controlled from PC

With the CS-2000/CS-2000A set to Remote Mode, the flow of communication commands and data for taking a measurement and reading the measured spectral irradiance, chromaticity, and measurement conditions is as shown in Figure 1 below.

- The instrument should already be set to Remote Mode by sending the command "RMTS, 1".
- The measuring button on the instrument should be disabled by sending the command "MSWE, 0".

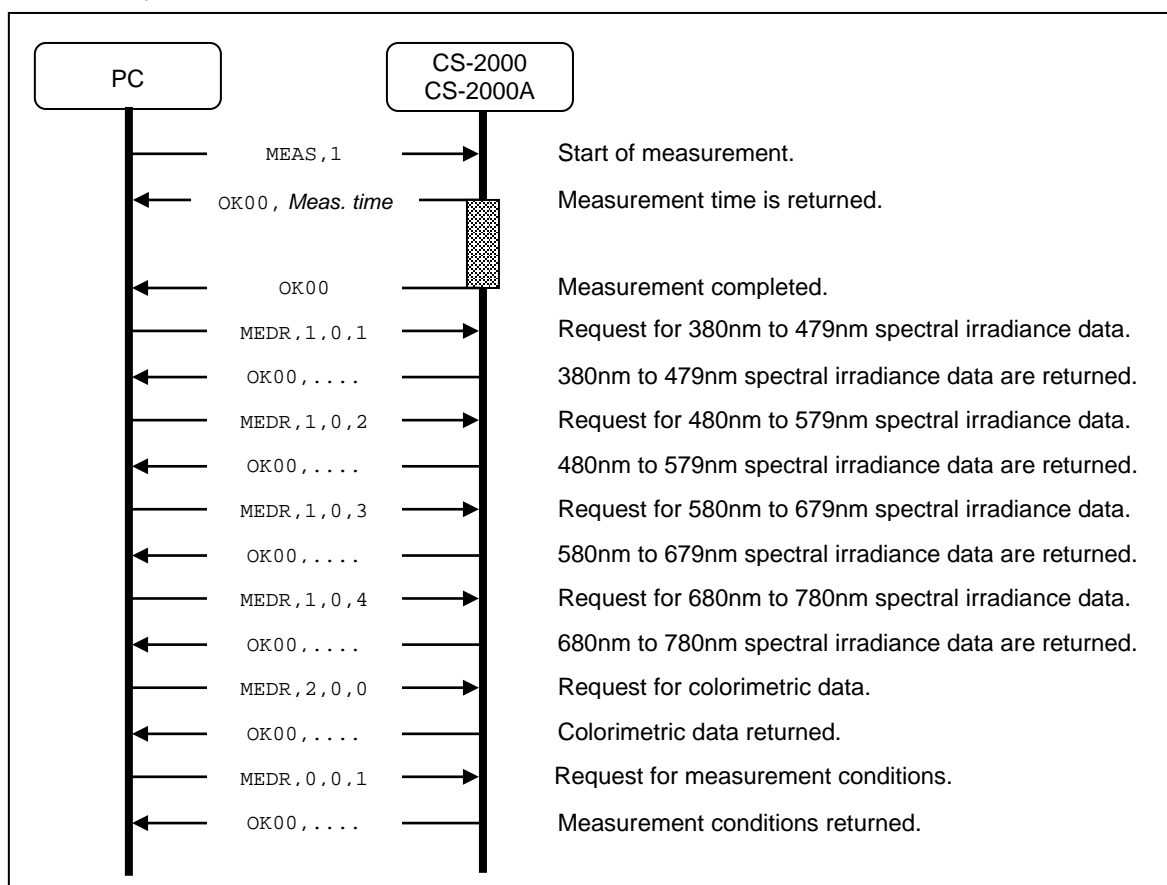


Figure 1: Communication flow from measurement to output of data and conditions



## 8.2 Communication flow for measurement using instrument measuring key

With the CS-2000/CS-2000A set to Remote Mode and the instrument measuring button enabled, the flow of operations and communication commands/data for taking a measurement and reading the measured spectral irradiance is as shown in Figure 2 below.

- The instrument should already be set to Remote Mode by sending the command "RMTS, 1".
- The measuring button on the instrument should be enabled by sending the command "MSWE, 1".
- When the instrument's measuring button is enabled, the data stored in the instrument buffer is cleared when either of the following actions is performed:
  - Reading of the spectral irradiance data for all wavelengths is completed.
  - Reading of any set of colorimetric data is completed.

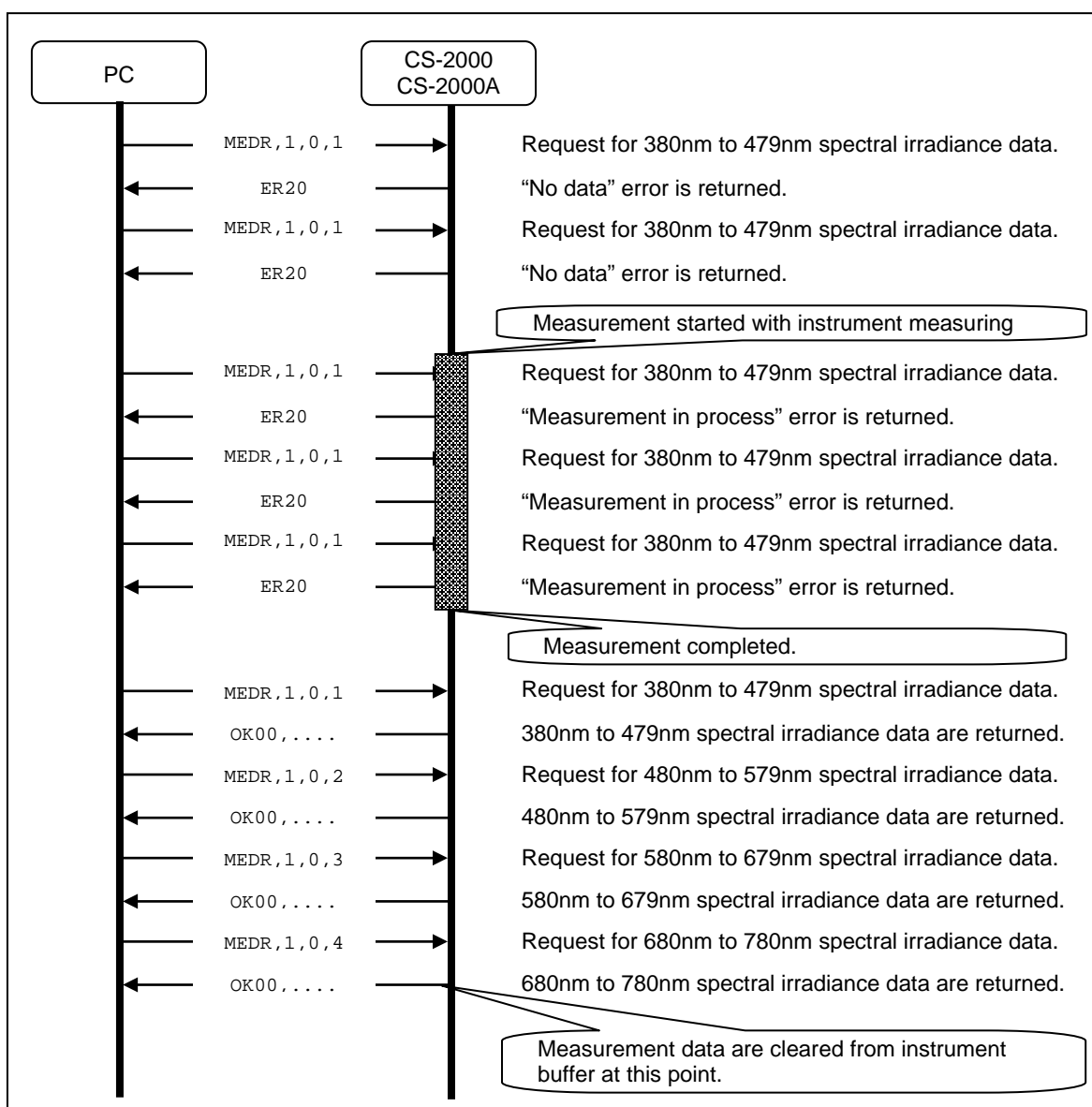


Figure 2: Communication flow for measurement using instrument measuring button

## 9 Modification of programs written for CS-2000 units with previous firmware versions

Modifications to programs written for CS-2000 units with firmware ver. 1.01.0000 or earlier in order to use the program with CS-2000 units with firmware ver. 1.10.0003 or later, or with CS-2000A units are described in the following sections.

### 9.1 Function changes in updated firmware

There are some function additions and changes in CS-2000/CS-2000A firmware ver. 1.10.0003 or later compared to ver. 1.01.0000 or earlier.

Function changed/added	Explanation
Speed Mode Read/Set ("SPMR"/"SPMS"): <ul style="list-style-type: none"> <li>Mode name change</li> <li>Additional mode</li> </ul>	<ul style="list-style-type: none"> <li>Name of previous "Multi Integ" mode changed to "Multi Integ Normal"</li> <li>"Multi Integ Fast." mode added</li> </ul>
Speed Mode Read/Set ("SPMR"/"SPMS"): <ul style="list-style-type: none"> <li>Function added to allow selection of "Internal ND Filter Mode" of "Off", "On", or "Auto" for "Speed Mode" of other than "Manual".</li> </ul>	When "Speed Mode" is set to "Normal", "Fast", "Multi Integ Normal" or "Multi Integ Fast", "Internal ND Filter Mode" can be set to "Off", "On", or "Auto".
Change of output for calculation error when reading measurement data, stored data, or target data.	When using communication commands to read out measurement data, the output data for a calculation error is the calculation error number as described in the Numerical Output Formats section. (For firmware ver. 1.10.0003 or earlier, "0" was output.)
Change in Identification Data Read ("IDDR")	Addition of a value for "Variation code" to indicate CS-2000A.

### 9.2 Command changes to correspond to updated firmware

Some communication commands and/or their parameters have been changed or added to correspond to the changes in function in the CS-2000/CS-2000A firmware ver. 1.10.0003 or later compared to ver. 1.01.0000 or earlier.

#### 9.2.1 Command changes to correspond to change of output for calculation error number

Function	Command	Details of change
Measurement mode setting	SPMS	When Command Parameter [1] is set to "0" (Normal) or "1" (Fast): Command Parameter [2] (Internal ND filter mode) must be added.
		When Command Parameter [1] is set to "2"(Multi Integ Normal) or "4"(Multi Integ Fast): Command Parameter [3] (Internal ND filter mode) must be added.
	SPMR	When Response Parameter [1] is set to "0" (Normal) or "1" (Fast): Output of Response Parameter [2] (Internal ND filter mode) has been added. Processing to handle the output of this added parameter should be added if necessary.
		When Response Parameter [1] is set to "2"(Multi Integ Normal) or "4"(Multi Integ Fast): Output of Response Parameter [3] (Internal ND filter mode) has

		been added. Processing to handle the output of this added parameter should be added if necessary.
	MEDR STDR TGDR	Upper limit of Response Parameter [1] is changed to "4" (Multi Integ Fast". If a check for the upper limit of Response Parameter [1] is currently implemented in the program, the upper limit should be changed accordingly.

### 9.2.2 Command changes to correspond to calculation error number changes

Function	Command	Details of change
Measurement Data Read Stored Data Read Target Data Read	MEDR STDR TGDR	If a calculation error occurs, the calculation error number will be output. Processing to handle the output of this calculation error number should be added. <ul style="list-style-type: none"> <li>For further information, refer to the Numerical Output Formats section.</li> </ul>

### 9.2.3 Command changes to correspond to Identification Data changes

Function	Command	Details of change
Instrument ID data Read	IDDR	A value of "2" (CS-2000A) may be output for Response Parameter [2]. Processing to handle the output of this value should be added if necessary.